

How to change a Marvel hearing aid battery with a tamperproof lock



1. Equipment you will need



2. There is a tiny hole normally at the bottom of the hearing aid to put the unlocking tool into.



3. Open up the battery door. Using the opening tool as a lever.



4. You will see the battery



5. Take out the dead battery.


You may need to push it out with a pencil or pen.



6. Put the dead battery back in pack and store out of reach of children. **Do not** put the sticker on it.



7. Batteries may be in an orange, brown or blue packet depending on type of hearing aid

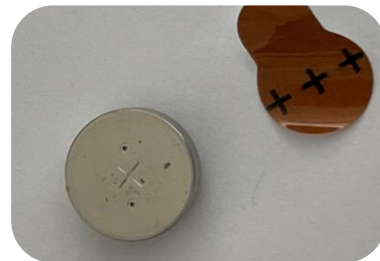

8. Make sure you have the right type of battery. Take the sticker off and wait for a few seconds.



9. Insert battery with flat side up so you can see the cross. There is usually a picture in the battery compartment to help.



10. Close the battery door. If it is difficult to close the battery door the battery may be in the wrong way.

Do not force the door as this can damage the hearing aid.



11. Check hearing aid is working in cupped hand. (you will hear a whistling noise)