



DAILY HEARING AID CHECKS



1. Turn the hearing aid on

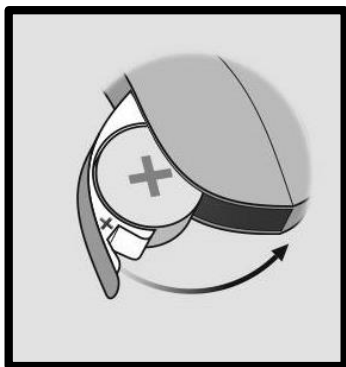

2. Cup it in the hands as shown

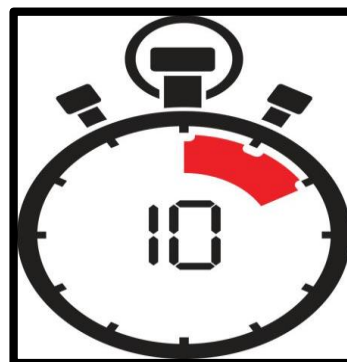

3. Count to ten to allow the aid to turn on and power up



4. The hearing aid should begin to whistle




5. Once it starts to whistle count to ten again.

If it keeps whistling, the battery is fine and the hearing aid is working; it is fine to use.




6. If there is no whistling sound, cup the hands more closely around the hearing aid.

If there is a sound, go back to step 4.

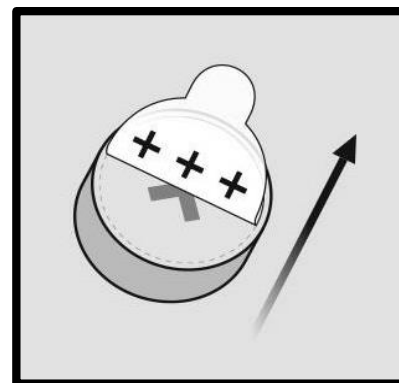

7. If there is still no sound, change the battery.

If there is a sound, go back to step 4.

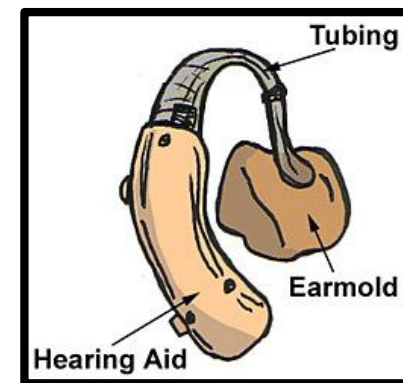

8. Check if the tubing or ear mould is blocked with wax – remove if possible.

Still no sound? Contact the child's carers and/or ETHV staff for advice

MY DAILY EQUIPMENT CHECKS

NAME:

WEEK BEGINNING:



CHECK	MONDAY	TUESDAY	WEDNESDAY	THURSDAY	FRIDAY
Batteries (working?) (batteries changed?)					
Tubing Check – retube if hard, discoloured or blocked.					
Moulds Cleaned with a wipe or in warm soapy water?					
Puffer Puff out water/bubbles from inside the tubing.					
Which Ear? Have I put my aid back in the right ear? Red = right Blue = left					
LING SOUND check Note down any difficulties or changes					
Additional Comments					