


# Cleaning Guide

## Slim tube

### How to clean your hearing aid slim tube




### Cleaning the dome

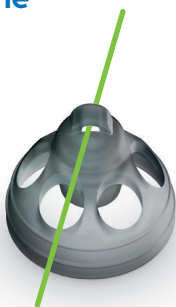

1. Push the green cleaning rod horizontally through the wax guard of the dome.
2. Clean the green cleaning rod with a damp cloth before removing.

## Cleaning the hook

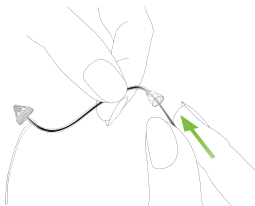

1. Unscrew the hook from the hearing aid anticlockwise.
2. Thread the green cleaning rod from the top to the bottom of the hook, stopping as it reaches the wax guard of the dome.
3. Clean any debris from the wax guard as described previously.
4. To re-attach the slim tube, screw back onto the hearing aid clockwise.

## Contact details

03033 306649

[nlg-tr.audiology@nhs.net](mailto:nlg-tr.audiology@nhs.net)

[www.nlg.nhs.uk/services/audiology/](http://www.nlg.nhs.uk/services/audiology/)